


Flag-Tag Mouse Monoclonal Antibody

Cat. No: DTF01

DESCRIPTION

Product Name	Flag-Tag Mouse Monoclonal Antibody
Host Species	Mouse
Purified by	Protein A affinity column
Applications	WB/IP
clone	M2, monoclonal
Specificity	DYKDDDDK Tag Antibody, mAb, Mouse recognizes C-terminal, N-terminal and internal tagged fusion proteins.
Immunogen	DYKDDDDK synthetic peptide conjugate to KLH
Ig class	IgG1
Conjugate	Unconjugated
Concentration	1.0mg/ml
Storage	For up to 12 months at -80℃. After receive, the antibody should be divided into the right amount. Avoid repeated freezing and thawing cycles.
Storage Buffer	PBS, pH 7.4
Recommended Dilutions	WB: 1 : 500-1 : 10000 IP: 1 : 200-1 : 2000

APPLICATION

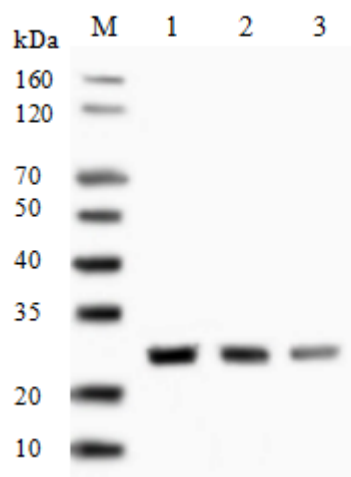

Lane 1: Flag Tagged Protein (0.80 μ g)

Lane 2: Flag Tagged Protein (0.50 μ g)

Lane 3: Flag Tagged Protein (0.10 μ g)

M: Western Blot Marker

Using Anti- Flag antibody