


His-Tag Mouse Monoclonal Antibody

Cat. No: DTH01

DESCRIPTION

Product Name	His-Tag Mouse Monoclonal Antibody
Host Species	Mouse
Purification	>95% Protein A purified
Applications	WB / ELISA / IF
Specificity	His-Tag Mouse Monoclonal antibody detects endogenous levels of internal, C-terminal, or N-terminal 6His-tagged proteins.
Ig class	Immunogen-specific mouse IgG1
Content	1.0 mg/ml with PBS, pH 7.4
Storage	For up to 12 months at -80°C. After receive, the antibody should be divided into the right amount. Avoid repeated freezing and thawing cycles.
Recommended Dilutions	WB: 1 : 2500-1 : 5000 ELISA: 1 : 2500-1 : 10000 IF: 1 : 500

APPLICATION

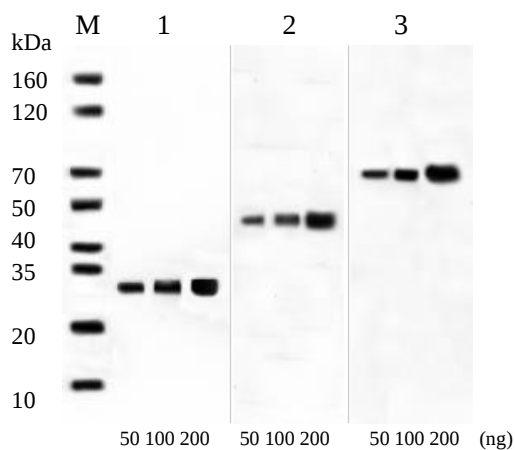

Western Blot analysis of HHHHHH tagged-fusion protein

Lane 1: C-His fusion protein


Lane 2: N-His fusion protein

Lane 3: Internal-His fusion protein

M: Molecular weight markers

Using His-Tag Mouse Monoclonal Antibody

(Dilute Proportion: 1 : 3000)

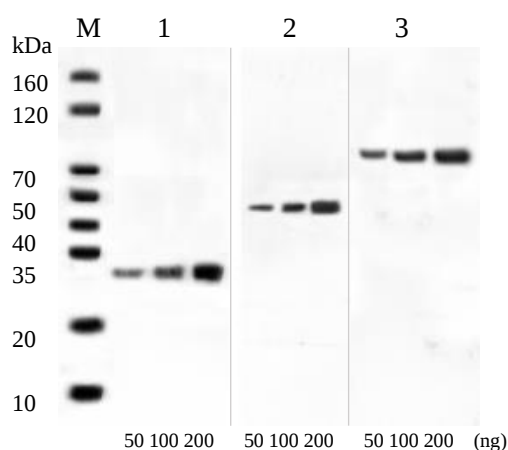

Western Blot analysis of HHHHHH tagged-fusion protein

Lane 1: C-His fusion protein


Lane 2: N-His fusion protein

Lane 3: Internal-His fusion protein

M: Molecular weight markers

Using His-Tag Mouse Monoclonal Antibody

(Dilute Proportion: 1 : 5000)

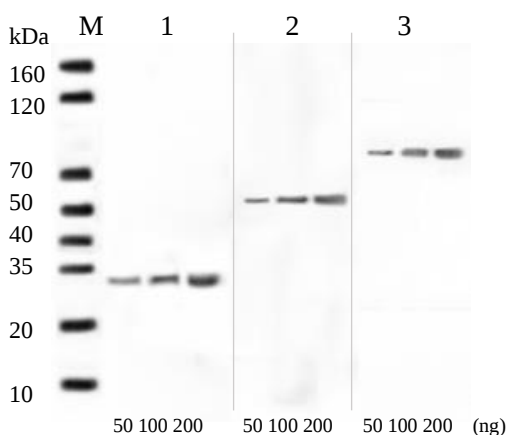

Western Blot analysis of HHHHHH tagged-fusion protein

Lane 1: C-His fusion protein

Lane 2: N-His fusion protein

Lane 3: Internal-His fusion protein

M: Molecular weight markers

Using His-Tag Mouse Monoclonal Antibody

(Dilute Proportion: 1 : 15000)