

Rabbit Anti TAP Tag Antibody, pAb

Cat.No: DTA0190

DESCRIPTION

Name	Rabbit Anti TAP Tag Antibody, pAb
Immunogen	TAP tag peptide conjugated to keyhole limpet haemocyanin(KLH)
Host Species	Rabbit
Specificity	This antibody recognizes TAP tag fusion proteins.
Form	Liquid
Concentration	0.5 mg/ml
Storage Buffer	1×PBS, pH7.4, containing 50% glycerol and 0.03% Proclin 300.
Purification	Antigen affinity purification
Clonality	Polyclonal
Isotype	Rabbit IgG
Shipping and Storage	Shipped in ice bag. The antibody can be stored for 1-2 weeks at 2-8°C. For long term storage, store it at -20°C. Avoid repeated freezing and thawing cycles.
Expiration Date	Specified on product label.

Applications

Application	Recommended Usage
ELISA	0.005-0.5μg/ml
Western Blot	0.1-1μg/ml
Other applications	User-optimized

Examples

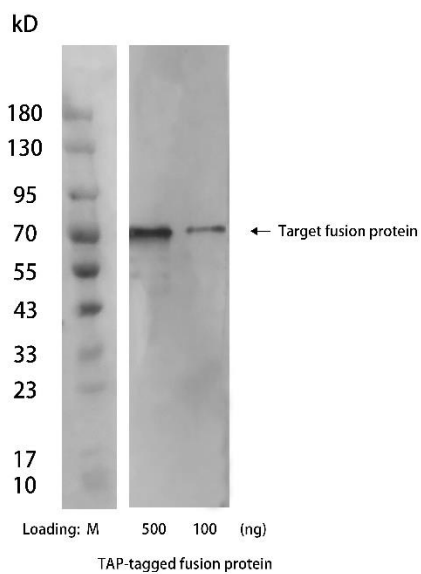


Figure 1. Detection of TAP tag in recombinant fusion protein were analyzed by Western blot using Rabbit Anti TAP Tag Antibody, pAb (Cat.DTA0190, at 1 μ g/mL concentration).

The signal was developed with Goat Anti-Rabbit IgG(H&L)[HRP], pAb(Cat.DTA0252).

感谢您选择德泰生物，我们竭诚为您服务！

For Research Use Only!