

Rabbit Anti DYKDDDDK Tag Antibody, pAb

Cat.No: DTF0201

DESCRIPTION

Name	Rabbit Anti DYKDDDDK Tag Antibody, pAb
Immunogen	DYKDDDDK peptide conjugated to keyhole limpet haemocyanin(KLH)
Host Species	Rabbit
Specificity	This antibody recognizes C-terminal, N-terminal, and internal region of DYKDDDDK tag fusion proteins.
Form	Liquid
Concentration	0.5 mg/ml
Storage Buffer	1×PBS, pH7.4, containing 50% glycerol and 0.03% Proclin 300.
Purification	Antigen affinity purification
Clonality	Polyclonal
Isotype	Rabbit IgG
Shipping and Storage	Shipped in ice bag. The antibody can be stored for 1-2 weeks at 2-8°C. For long term storage, store it at -20°C. Avoid repeated freezing and thawing cycles.
Expiration Date	Specified on product label.

Applications

Application	Recommended Usage
ELISA	0.025-0.1µg/ml
Western Blot	0.1-1µg/ml
IP	1-2µg/ml
Other applications	User-optimized

Examples

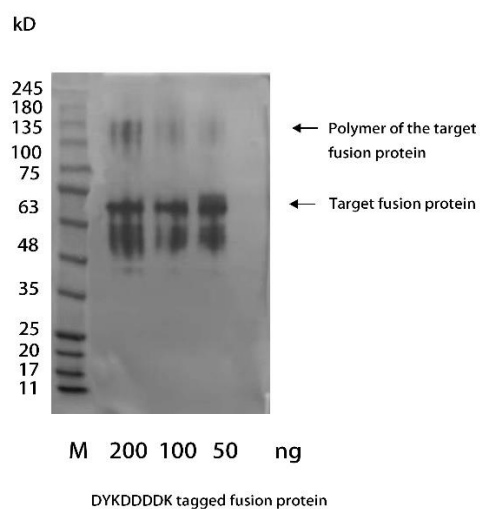


Figure 1. Detection of Flag tag in recombinant fusion protein were analyzed by Western blot using Rabbit Anti DYKDDDDK Tag Antibody, pAb(Cat.DTF0201, at 1 μ g/mL concentration).

The signal was developed with Goat Anti-Rabbit IgG(H&L)[HRP], pAb(Cat.DTA0252).

感谢您选择德泰生物，我们竭诚为您服务！

For Research Use Only!